


Pairing your Android Phone with your PC for usage with inSSIDer (Windows 7 Guide)

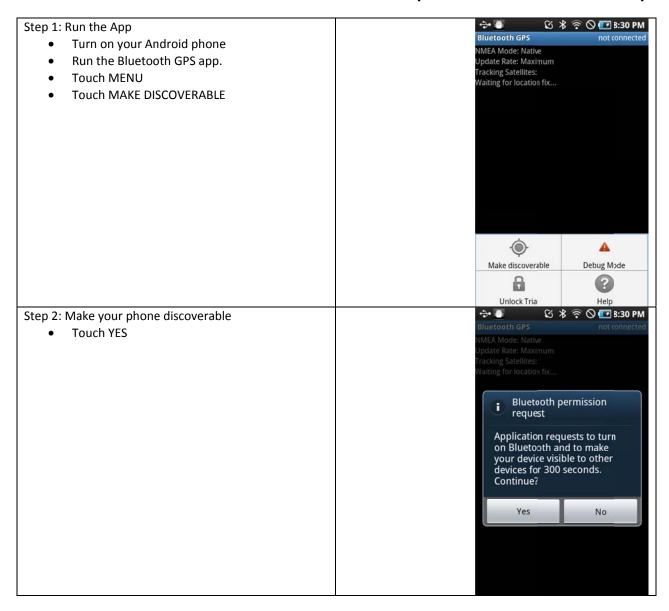

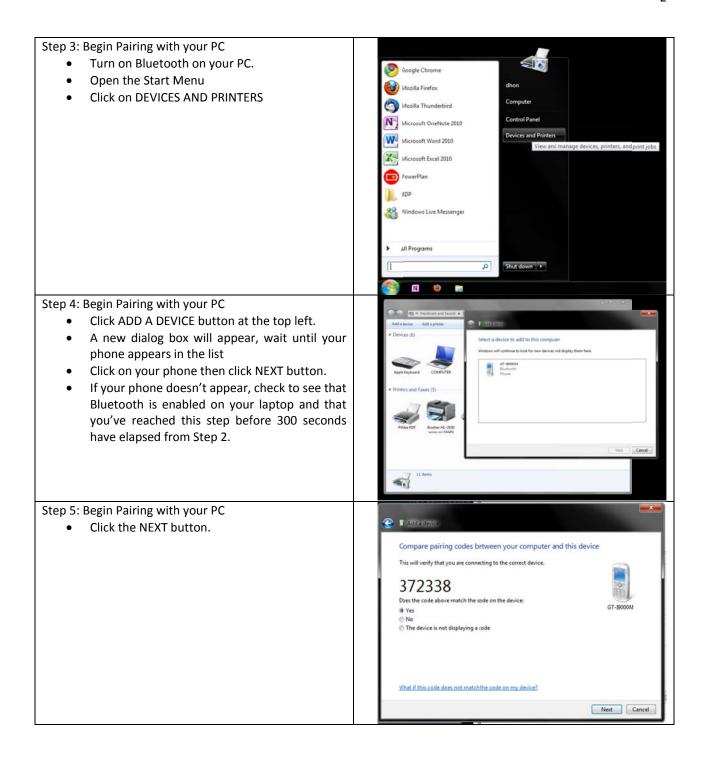

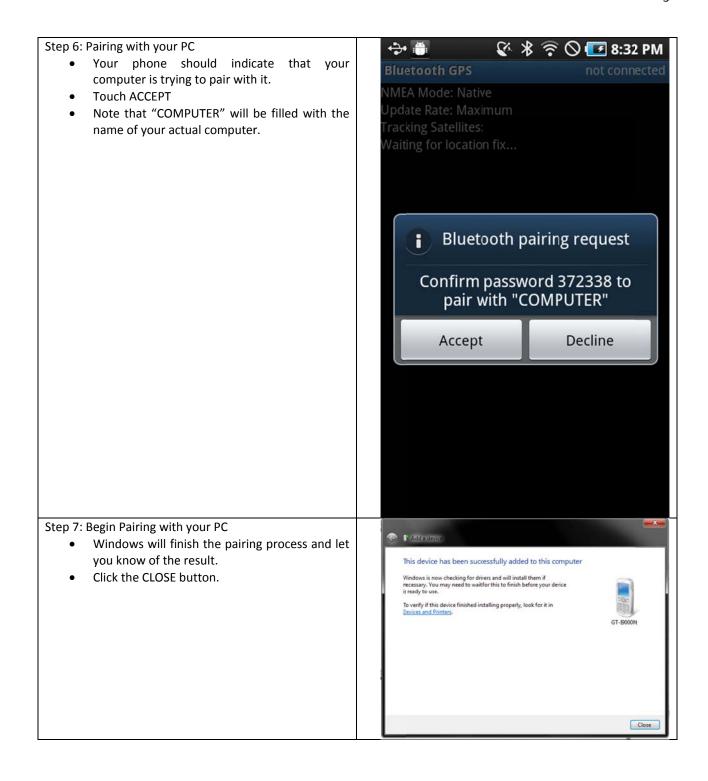

Step 8: Confirming the COM port number Look for your phone in the Device and Printer's Right click to display the context menu Left click on PROPERTIES Step 9: Confirming the COM port number GT-I9000M Properties The properties window will appear for your General Hardware Services Bluetooth phone. This Bluetooth device offers the following services. To use a service, select the check box. Click the SERVICES tab Look for "Serial port (SPP)" 'Android_Phone' or Bluetooth Services 'SPP Slave' COM4 ☑ Dial-up networking (DUN) Note the COM number to its right. Serial port (SPP) 'Android_Phone' COM5 In this example, COM5 is shown. Cancel Step 10: Configuring inSSIDer inSSIDer 2.0 Download and install inSSIDer Run inSSIDer File View Help Click FILE Configure GPS Click CONFIGURE GPS Change log filename Start Logging Convert GPX to KML Export to NS1 ondo Exit

Step 11: Configuring inSSIDer

- Pick the matching COM you noted from earlier in the PORT NAME setting.
- In this example, COM5 is the matching one for the Android Phone
- Note that if the COM number you are looking for does not appear in the list, simply type it in exactly. Example: COM5
- Click SAVE

Step 12: Configuring inSSIDer

- Click the STARTGPS at the top right side of inSSIDer window.
- Inssider will connect with your phone and pair match the detected WIFI networks and their approximate locations.

